

MALLARD SMALL ANIMAL COMBINED ANESTHESIA + VENTILATOR

[CLICK HERE TO REQUEST A QUOTE](#)


MODEL: 2200B

- Compact design suitable for medical facilities
- Includes the M2440V anesthesia ventilator
- 16 x 26 heavy shelf holds 2400V (see specs below) and absorber circuit
- Equipped with a small animal CO2 Absorber with a pressure gauge.
- Includes 19mm enclosed upright pop-off valve utilizing a proven unidirectional valve design
- Quick disconnect common gas outlet for convenient setup
- 0-4 LPM O2 flowmeter with mechanical stop flow control for precise oxygen delivery
- Oxygen flush function with safety restrictor ensures patient safety during emergencies
- Vaporizer mount in single or dual
- Quick soda lime canister change
- Absorber volume 2.5 liters
- Higher inspiratory flow rates for larger animals (up to 400 lbs)
- MRI Conditional upon request



2200B STANDARD SPECIFICATIONS

- **Classification:** A microprocessor based, electronic, time cycled, pressure limited, volume ventilator (MRI Optional)
- **Bellows Orientation:** Inverted (Ascending)
- **Volume:**
 - **Selectable** 200-2,200ml
 - **Optional:** Pediatric Bellows 50-300ml
 - **Optional:** Expanded Bellows 200-3,000ml
- **Controls (MRI Optional):**
 - **Inspiratory Time** - .1 - 3 seconds
 - **Inspiratory Flow** - 1 - 120 LPM
 - **Inspiratory Pressure** - 10 - 80cm H2O
 - **Rate** - 2 - 80 BPM
 - **Pneumatic Manual Push Button for Inspiration:** YES
- **LED Display (MRI Optional):**
 - Inspiratory Time
 - Rate (BPM)
 - I:E Ratio
- **Alarms and Alerts (MRI Optional):**
 - Failure to Cycle - Visual-Audible
 - Power Failure - Visual - Audible
 - Low Gas Pressure - Visual - Audible
- **Safety System:** Overpressure Governor 80cm H2O (adjustable)
- **Gas Scavenger:** 19mm port
- **Size:** 30" (76 cm) Width x 30" (76cm) Depth x 45"-61" (114-155 cm) Height
- **Weight:** 50 lbs (23 kg)
- **Options:** MRI Conditional

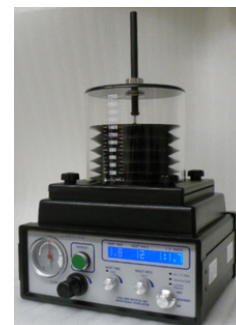

2400V



2400V-MRI

WARRANTY

Complete two year warranty

Mallard Medical

Phone: (530) 605-2522

Email: info@mallardmedical.com

Website: www.mallardmedical.com