

MALLARD SMALL ANIMAL ANESTHESIA VENTILATOR MRI CONDITIONAL

[CLICK HERE TO REQUEST A QUOTE](#)

MODEL: 2400V-MRI


- Ventilator has 100% pneumatic controls & displays no electronics in order to provide immunity to electronic emissions and interference.
- Contain all non-ferritic materials to avoid magnetic attraction.
- Inverted 2.2 liter bellows, ascending during exhalation as used in human anesthesia ventilators
- Patient may breathe spontaneously from the bellows
- Designed to adapt to small anesthesia machines
- Pneumatic display of selected respiratory variables
- Extremely quiet & smooth in operation

STANDARD SPECIFICATIONS

- **Classification:** A pneumatically controlled, time cycled, pressure limited, volume ventilator
- **Bellows Orientation:** Inverted (Ascending)
- **Volume:**
 - **Selectable** 200-2,200ml
 - **Optional: Pediatric Bellows** 50-300ml
 - **Optional: Expanded Bellows** 200-3,000ml
- **Controls:**
 - **Inspiratory Time** - .5 - 3 seconds
 - **Inspiratory Flow** - 1 - 120 LPM
 - **Inspiratory Pressure** - 10 - 80cm H₂O
 - **Rate** - 2 - 16 BPM
 - **Pneumatic Manual Push Button for Inspiration:** YES
- **Safety System:** Overpressure Governor 80cm H₂O (adjustable)
- **Gas Scavenger:** 19mm port
- **Size:**
 - 12" Width (30cm width)
 - 10" Depth (25 cm depth)
 - 19" Height (48 cm height)
- **Weight:** 10lbs (4.5 kg)



IDENTIFICATION SPECIFICATIONS

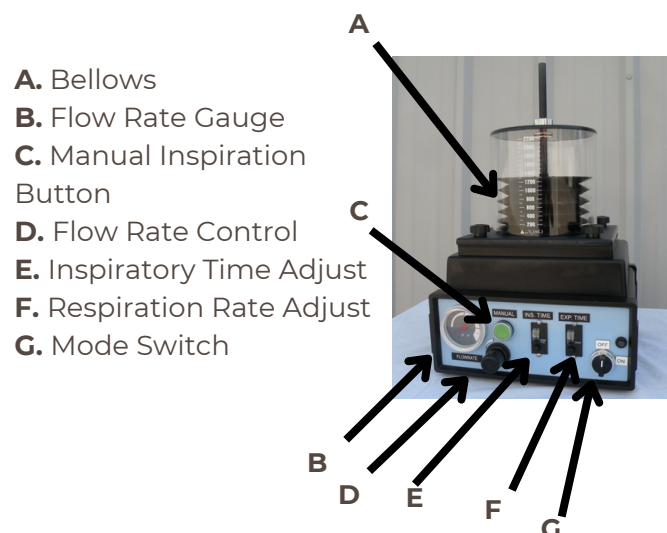

Mallard Medical

Phone: (530) 605-2522

Email: info@mallardmedical.com

Website: www.mallardmedical.com

WARRANTY

Complete two year warranty