

# MALLARD SMALL ANIMAL ANESTHESIA VENTILATOR

[CLICK HERE TO REQUEST A QUOTE](#)

## MODEL: 2400V


- Inverted bellows, ascending during exhalation as used in human anesthesia ventilators
- Patient may breathe spontaneously from the bellows
- Microprocessor based electronic control system
- L.C.D display of respiratory variables
- Flow rate adjustable from 0-120 liters per minute with expanded range of flow rate selectivity
- Designed to adapt to small animal anesthesia machines
- Extremely quite & smooth in operation.

## STANDARD SPECIFICATIONS

- **Classification:** A microprocessor based, electronic, time cycled, pressure limited, volume ventilator
- **Bellows Orientation:** Inverted (Ascending)
- **Volume:**
  - **Selectable** 200-2,200ml
  - **Optional: Pediatric Bellows** 50-300ml
  - **Optional: Expanded Bellows** 200-3,000ml
- **Controls:**
  - **Inspiratory Time** - .1 - 3 seconds
  - **Inspiratory Flow** - 1 - 120 LPM
  - **Inspiratory Pressure** - 10 - 80cm H2O
  - **Rate** - 2 - 80 BPM
  - **Pneumatic Manual Push Button for Inspiration:** YES
- **LED Display:**
  - Inspiratory Time
  - Rate (BPM)
  - I:E Ratio
- **Alarms and Alerts:**
  - Failure to Cycle - Visual-Audible
  - Power Failure - Visual - Audible
  - Low Gas Pressure - Visual - Audible
- **Safety System:** Overpressure Governor 80cm H2O (adjustable)
- **Gas Scavenger:** 19mm port
- **Size:** 12" Width x 10" Depth x 19" Height (30cm width) x (25 cm depth) x (48 cm height)
- **Weight:** 10lbs (4.5 kg)
- **Options:** MRI Conditional



## IDENTIFICATION SPECIFICATIONS

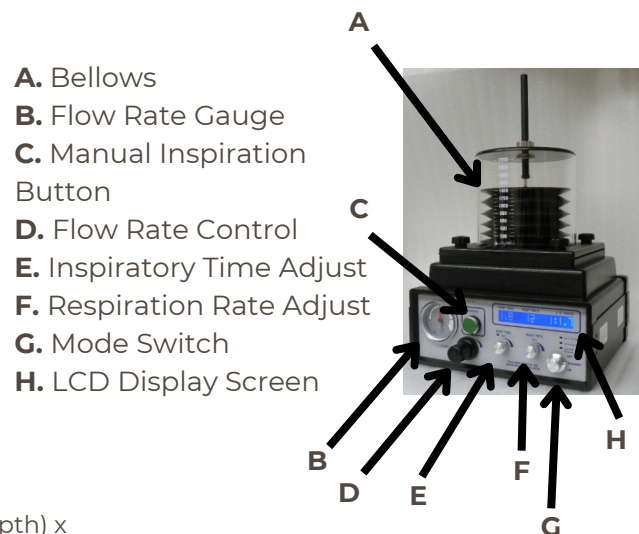

**Mallard Medical**

Phone: (530) 605-2522

Email: [info@mallardmedical.com](mailto:info@mallardmedical.com)

Website: [www.mallardmedical.com](http://www.mallardmedical.com)

### WARRANTY

Complete two year warranty